Sicherheitsinformationen

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Algemene veiligheidssymbolen:

Safety Instructions for the eaglefit HOME EMS System

EN

General Safety Instructions:

- Start your training at a low intensity and gradually increase the stimulation.
- Be sure to follow the instructions in the user manual.
- Do not use the device in rain, water, or damp environments.
- Children under 14 years old are not allowed to use the product.
- Do not use the product while operating a vehicle or machinery.
- Never wash the POWERBox in a washing machine and keep it away from water.

Mechanical Safety:

- Be aware of moving parts. Keep hair and clothing away to prevent injuries.

Electrical Safety:

- Never use the device in a wet environment.
- Ensure proper use to avoid electric shocks.
- Do not pinch or bend the cables and do not place heavy objects on them.
- Remove all piercings or similar items from the training area before training.
- If you find any damage, do not use the product and immediately inform the manufacturer.

Chemical Safety:

- Inquire about the materials used (silicone, nylon, elastane, carbon fibers) and be aware of possible allergies.

Thermal Safety:

- Never cover the POWERBox during operation to avoid overheating.


Instructions:

- Read the user manual carefully, pay attention to the contraindications, and watch the provided videos.

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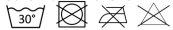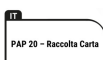
Page 1 of 8

MSDS Report

Product Name: EMS Machine (Containing Li-Ion Battery)

Model: HOME EMS, PRO EMS, STUDIO EMS, PANTS EMS, BELT EMS

Applicant: eaglefit GmbH

Written by: Sam Liu

Inspected by: Yemig

Approved by: Ray Zhou

Prepared By

East Notice Certification Service Co., Ltd.

1/F, Haohui Commercial Building, Zhuji Street, Dongpu Town, Tianhe District, Guangzhou City, China

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be accessible and confirmed at <http://www.enc-lab.com>.

East Notice Certification

1/F, Haohui Commercial Building, Zhuji Street,
Dongpu Town, Tianhe District, Guangzhou City

Tel: +86-020-8256 8534
E-mail: enc@enc-lab.com

Fax: +86-020-8256 8534
[Http:// www.enc-lab.com](http://www.enc-lab.com)

MSDS Report

Date: Mar. 31, 2023

Report No.: ENC2303232GZ94M

Page 2 of 8

MATERIAL SAFETY DATA SHEET (MSDS)

SECTION 1 – PRODUCT AND MANUFACTURER INFORMATION

1.1 Product identifiers

Product Name: EMS Machine (Containing Li-ion Battery)
Model: HOME EMS, PRO EMS, STUDIO EMS, PANTS EMS, BELT EMS
Brand: eaglefit

1.2 Relevant identified uses of the substance or mixture and uses advised against

Identified uses: Fitness and slimming

1.3 Details of the supplier of the safety data sheet

eaglefit GmbH
Add.: Tuchmacherstrasse 9, 89129 Langenau, Germany
Emergency telephone number: Tuchmacherstrasse 9, 89129 Langenau, Germany
Fax: 0049734591880065

SECTION 2 – HAZARDS IDENTIFICATION

2.1 Classification of the substance or mixture

This mixture does not meet the classification criteria according to OSHA.
This product is considered safe under normal use conditions.
The batteries are contained in the product. Thus, the ingredients have no hazard potential, except the battery is violated or dismantled.
If in case of mistreatment the ingredients are released, a spontaneously flammable gas mixture may be released under certain circumstances (measures according to SECTION 4 to 6).

2.2 GHS label elements, including precautionary statements

Pictogram(s)	No symbol.
Signal word	No signal word.
Hazard statement(s)	none
Precautionary statement(s) Prevention	none
Response	none
Storage	none
Disposal	none

2.3 Other hazards

none

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be accessible and confirmed at <http://www.enc-lab.com>.

ENC
East Notice Certification

1/F, Haohui Commercial Building, Zhuji Street,
Dongpu Town, Tianhe District, Guangzhou City

Tel: +86-020-8256 8534
E-mail: enc@enc-lab.com

Fax: +86-020-8256 8534
[Http:// www.enc-lab.com](http://www.enc-lab.com)

MSDS Report

Date: Mar. 31, 2023

Report No.: ENC2303232GZ94M

Page 3 of 8

SECTION 3 – INFORMATION ON INGREDIENT

The Li-ion Battery is installed in the product. Thus, the ingredients have no hazard potential, except the battery is violated or dismantled.

Ingredients of the Li-ion Battery:

Material	CAS No.	Weight (%)
The enclosure of the product is made of metal, with built-in lithium-ion battery		
Ingredients of the Li-ion Battery:		
Lithium	7439-93-2	2~3
Manganese dioxide	1313-13-9	30~40
Graphite	7782-42-5	3~4
Lithium perchlorate	7791-03-9	1
Ethylene glycol dimethyl ether	110-71-4	5.3
Propylene	108-32-7	5~8
Carbonate	7439-89-6	25~30
Chromium	7440-47-3	8~10
Molybdenum	7439-98-7	0.8~1.2
Polypropylene	9003-07-0	2

SECTION 4 – FIRST AID MEASURES

4.1 Description of first aid measures

- If inhaled:** Fumes can cause respiratory irritation. Remove to fresh air and consult a physician.
- In case of skin contact:** Immediately flush skin with plenty of water. If itch or irritation by chemical burn persists, consult a physician.
- In case of eye contact:** Immediately flush eye with plenty of water for at least 15 minutes. Consult a physician immediately.
- If swallowed:** If swallowing a battery, consult a physician immediately.
If contents come into mouth, immediately rinse by plenty of water and consult a physician.

4.2 Most important symptoms and effects, both acute and delayed

Transfer the patient to a safe place. Consult a doctor. Show this chemical safety technical manual to the doctor on site.

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC; this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be accessible and confirmed at <http://www.enc-lab.com>.

ENC
East Notice Certification

1/F, Haohui Commercial Building, Zhuji Street,
Dongpu Town, Tianhe District, Guangzhou City

Tel: +86-020-8256 8534
E-mail: enc@enc-lab.com

Fax: +86-020-8256 8534
[Http:// www.enc-lab.com](http://www.enc-lab.com)

MSDS Report

Date: Mar. 31, 2023

Report No.: ENC2303232GZ94M

Page 4 of 8

4.3 Indication of any immediate medical attention and special treatment needed

No data available

SECTION 5 – FIRE FIGHTING MEASURES

5.1 Extinguishing media

Suitable extinguishing media:

Extinguisher of alkaline metal fire is effective.

Plenty of cold water is also effective to cool the surrounding area and control the spread fire. But hydrogen gas may be evolved by the reaction of water and lithium and it can form an explosive mixture. Therefore in the case that lots of Li-ion Battery are burning in a confined space, use a smothering agent (e.g. carbon dioxide or dry sand).

5.2 Special hazards arising from the substance or mixture

Li-ion Battery may explode and release metal parts during fire.

5.3 Advice for firefighters

Use self contained breathing apparatus and full protective gear not to inhale harmful gas.

5.4 Further information

No data available

SECTION 6 – ACCIDENTAL RELEASE MEASURES

6.1 Personal precautions, protective equipment and emergency procedures

Wear personal protective equipment adapted to the situation (protection gloves, cloth).

6.2 Environmental precautions

Dispose off according to the local law and rules. Avoid leached substances to penetrate into the earth, canalization or water.

6.3 Methods and materials for containment and cleaning up

If the battery releases liquid, wipe it with a dry cloth.

6.4 Reference to other sections

For disposal see SECTION 13.

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be accessible and confirmed at <http://www.enc-lab.com>.

1/F, Haohui Commercial Building, Zhuji Street,
Dongpu Town, Tianhe District, Guangzhou City

Tel: +86-020-8256 8534
E-mail: enc@enc-lab.com

Fax: +86-020-8256 8534
[Http:// www.enc-lab.com](http://www.enc-lab.com)

MSDS Report

Date: Mar. 31, 2023

Report No.: ENC2303232GZ94M

Page 5 of 8

SECTION 7 – HANDLING AND STORAGE

7.1 Precautions for safe handling

In case the products are contained in equipment it is strongly recommended to follow the instruction manual requirements:

The product is required to be marked with a graphical symbol that alerts the user to refer to the instruction manual.

The instruction manual itself is required to contain

- a warning marking with text to alert the user of the potential chemical burn hazard associated with lithium batteries ingestion,
- an instruction as to the presence of a Lithium cell battery,
- possible effects of battery ingestion,
- an instruction to keep batteries away from children,
- an advice to seek immediate medical attention if it suspected that batteries have either been swallowed or placed inside any part of the body.

7.2 Conditions for safe storage, including any incompatibilities

Storage: Storage preferably at room temperature (approx. 20°C). Avoid large temperature changes. Do not store close to heating devices. Avoid direct sunlight. At higher temperature the electrical performance may be reduced. Preferred storage at 50% of the nominal capacity

Storage of unpacked batteries can cause short circuit and heat generation.

It is recommended to consider the "Technical Rule for Hazardous Substances TRGS 510 - Storage of hazardous substances in nonstationary containers" and to handle rechargeable Lithium ion batteries according to storage category 11 ("combustible solids").

SECTION 8 – EXPOSURE CONTROLS, PERSONAL PROTECTION

8.1 Control parameters

Keep away from heat and open flame.

Store in a cool dry place.

8.2 Exposure controls

Appropriate engineering controls: None required if used in a well-ventilated area.

Personal protective equipment:

- | | |
|--------------------------------|---|
| Eye Protection: | Not required during normal operations. |
| Hand Protection: | Protective gloves. |
| Body Protection: | Aprons are recommended. |
| Respiratory protection: | Respiratory protection is not required. |

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC; this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be processed and confirmed at <http://www.enc-lab.com>.

East Notice Certification

1/F, Haohui Commercial Building, Zhuji Street,
Dongpu Town, Tianhe District, Guangzhou City

Tel: +86-020-8256 8534
E-mail: enc@enc-lab.com

Fax: +86-020-8256 8534
[Http:// www.enc-lab.com](http://www.enc-lab.com)

MSDS Report

Date: Mar. 31, 2023

Report No.: ENC2303232GZ94M

Page 6 of 8

SECTION 9 – PHYSICAL AND CHEMICAL PROPERTIES

9.1 Information on basic physical and chemical properties

Appearance characters:	Rectangular solid
Nominal voltage:	3.7 V
Rated Capacity (mAh)	5000mAh

SECTION 10 – STABILITY AND REACTIVITY

10.1 Reactivity:	When heated above 100°C the risk of rupture occurs.
10.2 Chemical stability:	Stable (Performance deterioration depends on circumstance.)
10.3 Possibility of hazardous reactions:	Will not occur
10.4 Conditions to avoid:	See SECTION 7
10.5 Incompatible materials:	Water
10.6 Hazardous decomposition products:	Hydrogen

SECTION 11 – TOXICOLOGICAL INFORMATION

This product does not elicit toxicological properties during routine handling and use.
If the cells are opened through misuse or damage, discard recall. Internal components of cell are irritants and sensitizers.

SECTION 12 – ECOLOGICAL INFORMATION

When properly used or disposed, this product do not present environmental hazard.

SECTION 13 – DISPOSAL CONSIDERATIONS

When the battery is worn out, dispose of it under the ordinance of each local government or the law issued by relating government. Disposal of the worn-out battery may be subjected to Collection and Recycling Regulation.

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be accessible and confirmed at <http://www.enc-lab.com>.

East Notice Certification

1/F, Haohui Commercial Building, Zhuji Street,
Dongpu Town, Tianhe District, Guangzhou City

Tel: +86-020-8256 8534
E-mail: enc@enc-lab.com

Fax: +86-020-8256 8534
[Http:// www.enc-lab.com](http://www.enc-lab.com)

MSDS Report

Date: Mar. 31, 2023

Report No.: ENC2303232GZ94M

Page 7 of 8

SECTION 14 – TRANSPORT INFORMATION

UN Number	UN 3480 or UN 3481	
Transport hazard class(es)	9	
Packing group	--	
Marine pollutant	No	
UN Proper shipping name	☐ Lithium ion cell/battery = UN3480 with Section II of PI965; ☐ Lithium ion cell/battery packed with equipment = UN3481 with Section II of PI966; ☒ Lithium ion cell/battery contained in equipment = UN3481 with Section II of PI967.	
ICAO/IATA	ICAO-TI and IATA-DGR Li-Ion Battery according to NEW PACKING INSTRUCTION 965-967 of IATA DGR	DGR 64 th
IMDG CODE	International Maritime Dangerous Goods Code under Special Provision 188	IMDG CODE (Amdt.40-20)
ADR	European Agreement concerning the International Carriage of Dangerous Goods by Road under Special Provision 188	ADR 2023
RID	Regulations concerning the International Carriage of Dangerous Goods by Rail under Special Provision 188	RID 2023
DOT	The Transportation of Lithium-Ion cells and batteries are governed by US DOT CFR49 Part 171-180 of the US Hazardous Materials Regulations (HMR). CFR49 part 173.185(c) and the Special Provisions contained in 171.102 provide information on exceptions and packaging.	DOT 49 CFR
May be shipped without being declared as Class 9 dangerous goods, when meet to the requirements above. The dangerous goods regulations require that each battery design be subject to tests contained in Section 38.3 of the UN Manual of Tests and Criteria prior to being offered for transport.		

SECTION 15 – REGULATORY INFORMATION

International Conventions

* Air - IATA (International Air Transport Association) Dangerous Goods Regulations (DGR) 64th.

* Air - ICAO (International Civil Aviation Organization) Technical Instructions for the safe transport of dangerous goods by air.

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be accessible and confirmed at <http://www.enc-lab.com>.

1/F, Haohui Commercial Building, Zhuji Street,
Dongpu Town, Tianhe District, Guangzhou City

Tel: +86-020-8256 8534
E-mail: enc@enc-lab.com

Fax: +86-020-8256 8534
[Http:// www.enc-lab.com](http://www.enc-lab.com)

MSDS Report

Date: Mar. 31, 2023

Report No.: ENC2303232GZ94M

Page 8 of 8

* Sea – IMDG (International Maritime Dangerous Goods) IMDG CODE (Amdt.40-20). Under IMDG CODE Special Provision 188.

* Land – ADR 2023 (road) , RID 2023 (rail)

United Nations *Recommendations on the Transport of Dangerous Goods, Model Regulations – Ref. ST/SG/AC.10/30/Rev.9

United Nations *Recommendations on the Transport of Dangerous Goods, Manual of Tests and Criteria – Ref. ST/SG/AC.10/1/Rev.22 (Vol.I):2021

USA

* Code of Federal Regulations (49CFR Ch. 1 & 173 -185)

Both IATA and ICAO Special Provision A88 and IMO Special Provision 188, are identical to the requirements

EU Directive 2006/66/EC

OSHA hazard communication standard (29 CFR 1910.1200)

SECTION 16 –OTHER INFORMATION

The information above is believed to be accurate and represents the best information currently available to us. However, we make no warranty of merchantability or any other warranty, express or implied, with respect to such information, and we assume no liability resulting from its use. Users should make their own investigations to determine the suitability of the information for their particular purposes. In no event shall the company be liable for any claims, losses, or damages of any third party or for lost profits or any special, indirect, incidental, consequential, or exemplary damages howsoever arising, even if the company has been advised of the possibility of such damages.

Production of MSDS proving UN Manual of Tests and Criteria, Part III, sub-section 38.3 is met.

No	Item	Result	Remark
1	Altitude simulation	PASS	
2	Thermal test	PASS	
3	Vibration	PASS	
4	Shock	PASS	
5	External short circuit	PASS	
6	Impact	PASS	
7	Overcharge	PASS	
8	Force discharge	PASS	For cell only

MSDS Creation Date: Mar. 31, 2023

-----END OF REPORT-----

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be accessible and confirmed at <http://www.enc-lab.com>.

1/F, Haohui Commercial Building, Zhuji Street,
Dongpu Town, Tianhe District, Guangzhou City

Tel: +86-020-8256 8534
E-mail: enc@enc-lab.com

Fax: +86-020-8256 8534
[Http:// www.enc-lab.com](http://www.enc-lab.com)